


Power Floss®

DENTAL WATER JET

- A Reservoir
- B Max Fill Line
- C Tube
- D Button
- E Pump Housing
- F Nozzle

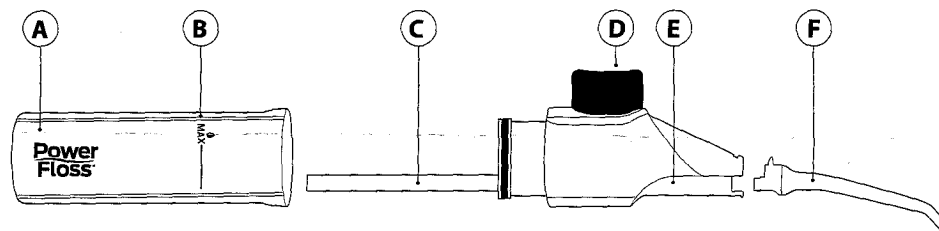

Prepare your Power Floss® Dental Water Jet

⚠ Important: Before using the Power Floss® Dental Water Jet, disconnect the Reservoir and submerge the Pump Housing, Tube, and Nozzle completely in water. Press the button for at least 60 seconds until it starts to bubble. Then lift the end of the Nozzle above the water and press the Button to ensure that the unit is working properly and has a steady jet stream.

Filling Your Power Floss® Dental Water Jet


1. Holding the Power Floss® Dental Water Jet upright, gently twist the Reservoir counter-clockwise and separate from the Pump Housing.
2. Fill the Reservoir to or below the max fill line on the Power Floss® Dental Water Jet body with either water or your choice of mouthwash.
Note: Only use mouthwash with the same consistency as water.
3. Reattach the Pump Housing by inserting it into the Reservoir and gently twist clockwise to lock into place.

Using Your Power Floss® Dental Water Jet

1. Position your thumb as shown in diagram 1 (see above), on the Button.
2. Always keep Power Floss® Dental Water Jet in an upright position during use to ensure a steady jet stream.
3. Press the Button several times until liquid comes out in a steady stream. You can adjust the speed and volume of the stream by varying the force you apply on the Button. For best results, experiment using different pressures and varying the angle of the Power Floss® Dental Water Jet as you move it across your mouth.

Diagram 1

⚠ When operating the Power Floss® Dental Water Jet DO NOT tilt beyond horizontal otherwise water will leak from the slits under the reservoir.

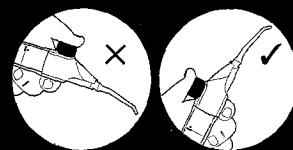

4. To clean in-between teeth: position the Nozzle against the teeth at the gum line, close your mouth around the Nozzle and press the Button. Use at both the front and back of the teeth for optimum results and continue for as long as required.
5. To clean bridges, braces, implants or other appliances: position the Nozzle under the bridge or appliance and press the Button. For best results, irrigate each site from more than one angle whenever possible.

Caring For Your Power Floss® Dental Water Jet

1. After each use, empty the Reservoir and depress the Button to remove any remaining fluid from within the unit.
2. For hygienic reasons do not share your Power Floss® Dental Water Jet.
3. Store your Power Floss® Dental Water Jet in a safe place and ensure the Nozzle cannot be damaged.
4. When traveling, store the Power Floss® Dental Water Jet Nozzle in the Reservoir.

Cleaning Your Power Floss® Dental Water Jet

After use, empty the reservoir and depress the button to remove any remaining fluid from the unit. Flush the Power Floss® Dental Water Jet thoroughly using warm water and wipe dry.

**For More Tips, Tricks, Instructional Videos, And To Purchase Additional Nozzles (Part F)
Please Visit: www.PowerFlossTips.Com**